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1. Specifications and Scope of Supply

Item no.:	PM-1580
Dimensions:	200 x 140 cm
Material:	Cotton/natural ticking, latex free
Colour:	Cream
ILN:	4250355315809
Storage:	-5 – 40°C protect from sunlight

Scope of supply:

- Turnsheets 200 x 140 cm with immobilisation straps
- Fitted Sheet (PM-1584) 200 x 90 cm with sliding surface 120 x 65 cm for width 90 - 100 cm
- Immobilisation straps 75 cm
- Non-slip mat PM-7010
- Manual

Made in Germany

max.
200 kg

Fitted Sheet

2. Use / Placement

The turn sheet is used with patients who are particularly sensitive to pain and with heavy patients or ones requiring special care. The concept does not comprise any obligatory care techniques or methods which are used irrespective of the patient but adapts to a patient's personal capabilities. Moving and laterally repositioning patients is made much easier. By rolling up the turn sheet against the patient's body in a lateral position and then fixing it in place with a bed sheet provides simple and gentle lateral position support.

The turn sheet consists of a large blanket and a fitted sheet, which is placed round the mattress and then fixed to the underside with fasteners. The upper side of the fitted sheet is coated with low-friction material. It is placed on the mattress in the patient's back area with the sliding surface face up. For protection this surface should be covered with an additional fitted bed sheet beforehand.


The underside of the blanket is also made of low-friction material whilst the upper side is made of breathable, highly skin-compatible, non-dyed natural ticking. There is a nonwoven between these two materials. This makes the turn sheet soft and adaptable. The blanket eventually lies on the sheet. The bed sheet and any incontinence liner can now be laid in place. Since the turn sheet is breathable, it can remain in the bed permanently. With the immobilisation straps the upper part can be fastened to

the top of the bed with the clips in order to prevent it from slipping down when the head section of the bed is raised. The turn sheet can also be used without the undersheet, depending on the patient's weight. By fitting the enclosed immobilisation mat (anti-slip mat) between the upper and lower sliding surfaces the slip is reduced or eliminated in this area.

When using the turn sheet with an alternating pressure mattress, please take note of the following:

As there are many different alternating pressure mattresses, the use of the turn sheet should be checked in advance.

This product is suitable for re-use.

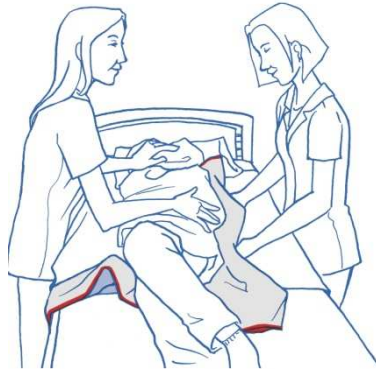3. Maintenance / Care

The product does not require any special maintenance. It can be washed in a washing machine at 90°C and is suitable for low dryer temperatures (45 – 50° C). When using non-standard washing agents, bleaching agents, softener or other chemical washing additives, the functional properties of the product can strongly degrade. Make sure that there is always enough ventilation between both sliding levels (danger of mould and humidity formation). The robust material can be disinfected with the usual wipe disinfectants. Check the compatibility before use where required on the basis of the disinfectant's instructions. We recommend the use of alcohol- and aldehyde-free agents, e.g. *Trionio®*. There is also some care information on the label.

4. Warranty

- Petermann gives a warranty that the product is devoid of material defects and manufacturing defects.
- The warranty is valid for a period of **3 years**. It does not cover any defects which are due to improper use.
- Petermann cannot be held liable for any accidental or indirect injury or damage.
- Make sure that there is always enough ventilation between both sliding levels (danger of mould and humidity formation).

5. Visual and Function Inspection / Safety Information / Re-use

Repositioning products are also subject to natural aging. Check the equipment regularly and carefully to make sure it is working properly and especially before re-use. Make sure that the product is in good condition and that the material and the seams are not damaged. Check the function of the product and make sure that the material resists the load. The smooth surface of the low-friction material makes it easy to slide and move. The best way to test these properties is by comparing the product with a new product.

For re-using the product it must be washed and disinfected before, see 3.

WARNING:

- ***Faulty products must no longer be used.***
- ***If the patient's mobility is restricted, it is essential to supervise the motions of transfer by providing trusted help.***
- ***The turn sheet may only be used in beds with protective side bars.***
- ***We recommend that the sheet only be used in combination with alternating pressure after prior testing. Contact our service team if necessary.***
- ***When using non-standard washing agents, bleaching agents, softener or other chemical washing additives, the functional properties of the product can strongly degrade.***
- ***The smooth surface considerably minimises resistance to friction during the change of position. Carefully check the transfer surfaces, making sure they are stable and brakes are applied (e. g. wheelchair, bed).***
- ***Make sure that there is always enough ventilation between both sliding levels (danger of mould and humidity formation).***

6. Product Labeling

The Petermann Alpha[®] turn sheet is a medical product class I and complies with the MDD directive 93/42/EEC.

Petermann GmbH is certified in accordance with the directive DIN EN ISO 13485 by the German TUV.

7. Disposal Information

For the turn sheet there are no special disposal rules; it can be disposed of in accordance with national and local regulations as necessary.

If you have any questions about the use of our products, please call our service number:

+49 (0) 9868 9339-0

Petermann GmbH • Ziegelhaus 12 • D-91601 Dombühl
Postal address: Petermann GmbH • Waldweg 3 • D-91601 Dombühl
Phone: +49 (0) 9868 9339-0 • Fax: +49 (0) 9868 9339-20
E-mail: petermann@pm-med.de • www.pm-med.de

"Information subject to technical modifications to our products and deviations in colour or printing errors"

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Petermann GmbH est certifié par le TUV conforme à la directive DIN EN ISO 13485.
Adresse postale: Petermann GmbH • Waldweg 3 • D-91601 Dombühl